

Premium Certified Reference Material RRM Fe-114

PRODUCT: IRON ORE Pilbara Region Western Australian CERTIFICATE OF ANALYSIS

Certificate Number: RRM CRM Fe-114 Rev 000

Date: 29 May 2017

Table 1 Fe-114 Constituents (ISO 2596:2006 Hygroscopic Moisture Corrected)

Analyte (Unit)	Assigned value	Standard Deviation		95% Confidence Limits		Number of Laboratories	Number of Analyses	Mean of number of analyses per Laboratory
		Within Laboratory	Between Laboratory	Lower	Upper			
Fe (%)	57,51	0,41	0,41	57,40	57,62	15	58	3,9
Fe calculated (%)	57,57	0,08	0,09	57,54	57,60	15	58	3,9
SiO ₂ (%)	5,17	0,05	0,06	5,16	5,19	15	58	3,9
Al ₂ O ₃ (%)	2,97	0,03	0,03	2,96	2,98	15	58	3,9
TiO ₂ (%)	0,102	0,002	0,004	0,101	0,104	15	58	3,9
Mn (%)	0,296	0,003	0,003	0,295	0,297	15	58	3,9
CaO (%)	0,052	0,008	0,009	0,049	0,055	15	58	3,9
MgO (%)	0,087	0,009	0,011	0,083	0,092	15	58	3,9
K ₂ O (%)	0,014	0,001	0,002	0,013	0,014	15	58	3,9
Na ₂ O (%)	0,028	0,007	0,012	0,022	0,034	13	51	3,9
P (%)	0,043	0,001	0,001	0,043	0,044	15	58	3,9
S (%)	0,033	0,001	0,002	0,032	0,034	14	54	3,9
LOI-371 %	7,52	0,06	0,13	7,42	7,62	5	30	6,0
LOI-425 %	7,77	0,04	0,05	7,74	7,81	8	44	5,5
LOI-650 %	8,34	0,04	0,07	8,30	8,38	11	54	4,9
LOI-1000 %	8,62	0,10	0,12	8,58	8,65	17	80	4,7
Cl (%)	0,002	0,001	0,002	0,001	0,002	12	48	4,0
As (%)	0,009	0,001	0,003	0,007	0,011	13	50	3,8
Ba (%)	-	-	-	-	-	11	45	4,1
Co (%)	0,002	0,001	0,001	0,001	0,002	12	49	4,1
Cu (%)	0,005	0,002	0,002	0,004	0,006	13	51	3,9
Cr (%)	0,005	-	0,001	0,004	0,005	11	36	3,3
Ni (%)	0,004	0,003	0,003	0,002	0,005	10	45	4,5
Pb (%)	0,003	0,001	0,001	0,002	0,000	10	43	4,3
Sn (%)	0,001	-	0,001	-	0,001	10	43	4,3
Sr (%)	0,001	0,001	0,001	-	0,001	9	39	4,3
V (%)	-	-	-	-	-	9	43	4,8
Zn (%)	0,001	0,001	0,001	0,001	0,002	10	42	4,2
Zr (%)	0,003	0,003	0,002	0,002	0,005	11	37	3,4

Table 2 Fe-114 Constituents (non-hygroscopically moisture corrected)

Analyte (Unit)	Assigned value	Standard Deviation		95% Confidence Limits		Number of Laboratories	Number of Analyses	Mean of number of analyses per Laboratory
		Within Laboratory	Between Laboratory	Lower	Upper			
Fe (%)	57,32	0,14	0,22	57,23	57,42	15	60	4,0
Fe calculated (%)	57,41	0,09	0,11	57,36	57,45	15	60	4,0
SiO ₂ (%)	5,15	0,03	0,04	5,14	5,16	15	60	4,0
Al ₂ O ₃ (%)	2,95	0,02	0,03	2,94	2,96	15	60	4,0
TiO ₂ (%)	0,101	0,001	0,002	0,100	0,102	15	60	4,0
Mn (%)	0,293	0,002	0,002	0,292	0,294	14	56	4,0
CaO (%)	0,049	0,003	0,006	0,046	0,052	15	60	4,0
MgO (%)	0,091	0,010	0,029	0,076	0,107	15	60	4,0
K ₂ O (%)	0,013	0,002	0,003	0,012	0,015	14	56	4,0
Na ₂ O (%)	0,020	0,009	0,012	0,014	0,025	11	42	3,8
P (%)	0,044	0,001	0,001	0,044	0,045	15	60	4,0
S (%)	0,032	0,001	0,003	0,030	0,034	14	56	4,0
LOI-371 %	7,76	0,03	0,15	7,66	7,87	9	45	5,0
LOI-425 %	8,06	0,07	0,10	8,00	8,12	7	39	5,6
LOI-650 %	8,58	0,06	0,08	8,54	8,62	11	56	5,1
LOI-1000 %	8,85	0,06	0,13	8,80	8,91	17	83	4,9
Cl (%)	0,002	-	0,001	0,002	0,003	9	36	4,0
As (%)	0,010	0,001	0,003	0,008	0,012	9	35	3,9
Ba (%)	0,002	0,001	0,002	0,001	0,004	6	26	4,3
Co (%)	0,003	0,001	0,002	0,002	0,004	10	37	3,7
Cu (%)	0,007	0,001	0,003	0,006	0,009	10	40	4,0
Cr (%)	0,007	0,001	0,003	0,005	0,009	8	30	3,8
Ni (%)	0,006	0,004	0,005	0,004	0,008	6	26	4,3
Pb (%)	0,004	0,001	0,001	0,004	0,005	6	23	3,8
Sn (%)	0,003	0,001	0,001	0,002	0,004	6	27	4,5
Sr (%)	-	0,001	0,001	-	-	5	18	3,6
V (%)	0,001	0,001	0,001	-	0,001	4	17	4,3
Zn (%)	0,002	0,001	0,001	0,001	0,003	6	23	3,8
Zr (%)	0,005	0,001	0,002	0,004	0,007	5	21	4,2

Introduction:

Certified Reference Materials (CRMs) are used by laboratories to prove the value of their service offerings and for clients of laboratories to evaluate and monitor laboratory performance. CRMs must comply with high metrological requirements and ensuring traceability of measurement results.

Since most techniques employing analytical instrumentation are comparative, these techniques require a sample of known composition (CRM) for accurate calibration. Grade and Matrix matched CRM's are thus vital to the core of the analytical chemistry industry.

Preparation of Material:

The material constituting RRM Fe-114 has been according to ISO 17034:2016, and includes the following:

- Drying to constant mass
- Crushing and dry milling to nominal 53µm particle size
- Homogenization
- Systematic rotary division of the entire lot, to final aliquot
- Packaging
- Rotary divided to nominal 250g sealed jars
- Rotary divided to nominal 10g geochem pouches, vacuum sealed in barrier foil
- Custom package sizes available on request

Unique to RRM, rotary division of all material to final packaging.

Methods of Analysis:

The analysis of the test samples has been conducted according to each individual laboratory's routine analytical procedures. The material has been dried at 105°C until constant mass was achieved, or corrected for moisture according to ISO 2596:2006. Each test sample has been prepared, analyzed and reported in duplicate; with results reported on a dry basis.

Elements and Oxides determined as follows:

- Multi-element Iron Ore Suite – XRF fused disc analysis
- Loss on Ignition – Thermo Gravimetric Analysis

Homogeneity Evaluation:

For the evaluation of the homogeneity of the Fe-114 material, 14 samples were selected throughout the batch for analysis in duplicate. ANOVA (one-way analysis of variance) was used to assess the homogeneity for measured Iron and total Loss on Ignition. No significant variation was observed for both 'batch' and 'within unit' exercises. The results and statistical evaluation from the interlaboratory studies (presented in Tables 3, 4, 5 and 6) further validate the homogeneity. The material can thus be considered fit for purpose for use as a CRM.

Statistical Evaluation:

We have grouped the laboratories according to the method of determining moisture. Certain laboratories determine the dry weight of iron ore samples at 105°C, whilst some use the ISO 2596:2006 correction for hygroscopic moisture. Some iron ore types report different moisture values for each of these methods, which in turn influences the chemical analysis as the determined moisture content has an influence on the dry starting mass.

A comprehensive statistical evaluation of the results received from the various laboratories was performed. Outliers were identified and removed from the data sets for the certification exercise.

The assigned values are the mean of means after removal of outliers. All constituents reported with a between laboratory % RSD of <5% are regarded as certified values.

Table 3 Fe-114 Statistical Evaluation (Hygroscopic Moisture Corrected)

Analyte (Unit)	Fe (%)	Fe calc (%)	SiO ₂ (%)	Al ₂ O ₃ (%)	TiO ₂ (%)	Mn (%)	CaO (%)
Assigned value	57.51	57.57	5.17	2.97	0.102	0.296	0.052
Number of Laboratories	15	15	15	15	15	14	15
Number of Analyses	58	58	58	58	58	56	58
Between Laboratory Standard Deviation	0.41	0.09	0.06	0.03	0.004	0.003	0.009
Within Laboratory Standard Deviation	0.41	0.08	0.05	0.03	0.002	0.003	0.008
Between Laboratory Relative Standard Deviation	0.71 %	0.16 %	1.08 %	1.07 %	3.45 %	1.02 %	16.9 %
Within Laboratory Relative Standard Deviation	0.71 %	0.13 %	1.01 %	1.00 %	1.81 %	0.87 %	14.9 %
Standard Uncertainty	0.055	0.016	0.009	0.005	0.001	0.001	0.002
Lower confidence limit	57.40	57.54	5.16	2.96	0.101	0.295	0.049
Upper confidence limit	57.62	57.60	5.19	2.98	0.104	0.297	0.055
Lower limit of tolerance	56.70	57.39	5.06	2.90	0.095	0.290	0.034
Upper limit of tolerance	58.33	57.75	5.28	3.03	0.109	0.302	0.069

Analyte (Unit)	P (%)	S (%)	MgO (%)	K ₂ O (%)	Na ₂ O (%)	As (%)	Cl (%)
Assigned value	0.043	0.033	0.087	0.014	0.028	0.002	0.009
Number of Laboratories	15	14	15	15	13	12	13
Number of Analyses	58	54	58	58	51	48	50
Between Laboratory Standard Deviation	0.001	0.002	0.011	0.002	0.012	0.002	0.003
Within Laboratory Standard Deviation	0.001	0.001	0.009	0.001	0.007	0.001	0.001
Between Laboratory Relative Standard Deviation	2.68 %	6.58 %	13.2 %	11.2 %	44%	102 %	36 %
Within Laboratory Relative Standard Deviation	1.98 %	4.12 %	9.81 %	8.08 %	25 %	49 %	15.1 %
Standard Uncertainty	-	0.001	0.002	-	0.003	-	0.001
Lower confidence limit	0.043	0.032	0.083	0.013	0.022	0.001	0.007
Upper confidence limit	0.044	0.034	0.092	0.014	0.034	0.002	0.011
Lower limit of tolerance	0.041	0.028	0.064	0.011	0.003	-	0.003
Upper limit of tolerance	0.046	0.037	0.110	0.017	0.052	0.005	0.016

Analyte (Unit)	Co (%)	Cu (%)	Pb (%)	Zn (%)	Ba (%)	Cr (%)	Ni (%)
Assigned value	-	0.002	0.005	0.005	0.004	0.003	0.001
Number of Laboratories	11	12	13	11	10	10	10
Number of Analyses	45	49	51	36	45	43	43
Between Laboratory Standard Deviation	-	0.001	0.002	0.001	0.003	0.001	0.001
Within Laboratory Standard Deviation	-	0.001	0.002	-	0.003	0.001	-
Between Laboratory Relative Standard Deviation	-	91 %	42 %	27 %	82 %	55 %	143 %
Within Laboratory Relative Standard Deviation	-	45 %	36 %	8.75 %	34 %	31 %	56 %
Standard Uncertainty	-	-	-	-	0.001	-	-
Lower confidence limit	-	0.001	0.004	0.004	0.002	0.002	-

Upper confidence limit	-	0.002	0.006	0.005	0.005	0.003	0.001
Lower limit of tolerance	-	-	0.001	0.002	-	-	-
Upper limit of tolerance	-	0.004	0.010	0.007	0.010	0.005	0.003

Analyte (Unit)	Sn (%)	Sr (%)	V (%)	Zr (%)
Assigned value	0.001	-	0.001	0.003
Number of Laboratories	9	9	10	11
Number of Analyses	39	43	42	37
Between Laboratory Standard Deviation	0.001	-	0.001	0.002
Within Laboratory Standard Deviation	0.001	-	0.001	0.001
Between Laboratory Relative Standard Deviation	193 %	-	56 %	74 %
Within Laboratory Relative Standard Deviation	182 %	-	42 %	18.9 %
Standard Uncertainty	-	-	0.000	0.001
Lower confidence limit	-	-	0.001	0.002
Upper confidence limit	0.001	-	0.002	0.005
Lower limit of tolerance	-	-	-	-
Upper limit of tolerance	0.003	-	0.003	0.008

Table 4 Fe-114 Statistical Evaluation (non-hygroscopically moisture corrected)

Analyte (Unit)	Fe (%)	Fe calc (%)	SiO ₂ (%)	Al ₂ O ₃ (%)	TiO ₂ (%)	Mn (%)	CaO (%)
Assigned value	57.32	57.41	5.15	2.95	0.101	0.293	0.049
Number of Laboratories	15	15	15	15	15	14	15
Number of Analyses	60	60	60	60	60	56	60
Between Laboratory Standard Deviation	0.22	0.11	0.04	0.03	0.002	0.002	0.006
Within Laboratory Standard Deviation	0.14	0.09	0.03	0.02	0.001	0.002	0.003
Between Laboratory Relative Standard Deviation	0.38 %	0.19 %	0.70 %	0.98 %	1.85 %	0.79 %	12.8 %
Within Laboratory Relative Standard Deviation	0.25 %	0.16 %	0.64 %	0.67 %	1.37 %	0.68 %	5.61 %
Standard Uncertainty	0.046	0.021	0.006	0.006	0.000	0.000	0.001
Lower confidence limit	57.23	57.36	5.14	2.94	0.100	0.292	0.046
Upper confidence limit	57.42	57.45	5.16	2.96	0.102	0.294	0.052
Lower limit of tolerance	56.89	57.18	5.08	2.89	0.097	0.288	0.036
Upper limit of tolerance	57.76	57.63	5.22	3.01	0.105	0.298	0.061

Analyte (Unit)	P (%)	S (%)	MgO (%)	K ₂ O (%)	Na ₂ O (%)	As (%)	Cl (%)
Assigned value	0.044	0.032	0.091	0.013	0.020	0.002	0.010
Number of Laboratories	15	14	15	14	11	9	9
Number of Analyses	60	56	60	56	42	36	35
Between Laboratory Standard Deviation	0.001	0.003	0.029	0.003	0.012	0.001	0.003

Within Laboratory Standard Deviation	0.001	0.001	0.010	0.002	0.009	-	0.001
Between Laboratory Relative Standard Deviation	2.51 %	9.94 %	32 %	22 %	60 %	41 %	33 %
Within Laboratory Relative Standard Deviation	1.19 %	2.27 %	10.8 %	15.6 %	45 %	21 %	7.06 %
Standard Uncertainty	-	0.001	0.008	0.001	0.003	-	0.001
Lower confidence limit	0.044	0.030	0.076	0.012	0.014	0.002	0.008
Upper confidence limit	0.045	0.034	0.107	0.015	0.025	0.003	0.012
Lower limit of tolerance	0.042	0.026	0.033	0.008	-	-	0.003
Upper limit of tolerance	0.046	0.038	0.149	0.019	0.043	0.004	0.017

Analyte (Unit)	Co (%)	Cu (%)	Pb (%)	Zn (%)	Ba (%)	Cr (%)	Ni (%)
Assigned value	0.002	0.003	0.007	0.007	0.006	0.004	0.003
Number of Laboratories	6	10	10	8	6	6	6
Number of Analyses	26	37	40	30	26	23	27
Between Laboratory Standard Deviation	0.002	0.002	0.003	0.003	0.005	0.001	0.001
Within Laboratory Standard Deviation	0.001	0.001	0.001	0.001	0.004	0.001	0.001
Between Laboratory Relative Standard Deviation	99 %	71 %	39 %	40 %	78 %	19.0 %	36 %
Within Laboratory Relative Standard Deviation	54 %	19.0 %	14.2 %	12.2 %	70 %	14.3 %	34 %
Standard Uncertainty	0.001	0.001	0.001	0.001	0.001	-	-
Lower confidence limit	0.001	0.002	0.006	0.005	0.004	0.004	0.002
Upper confidence limit	0.004	0.004	0.009	0.009	0.008	0.005	0.004
Lower limit of tolerance	-	-	0.002	0.001	-	0.003	0.001
Upper limit of tolerance	0.006	0.007	0.013	0.012	0.015	0.006	0.005

Analyte (Unit)	Sn (%)	Sr (%)	V (%)	Zr (%)
Assigned value	-	0.001	0.002	0.005
Number of Laboratories	5	4	6	5
Number of Analyses	18	17	23	21
Between Laboratory Standard Deviation	0.001	0.001	0.001	0.002
Within Laboratory Standard Deviation	0.001	0.001	0.001	0.001
Between Laboratory Relative Standard Deviation	148 %	92 %	62 %	37 %
Within Laboratory Relative Standard Deviation	148 %	92 %	39 %	17.7 %
Standard Uncertainty	-	-	-	0.001
Lower confidence limit	-	-	0.001	0.004
Upper confidence limit	-	0.001	0.003	0.007
Lower limit of tolerance	-	-	-	0.001
Upper limit of tolerance	0.002	0.002	0.004	0.009

LOI Certified Values:

The LOI statistics from the inter-laboratory study are presented in the table below:

Table 5 Fe-114 LOI Certified Values (Hygroscopic Moisture Corrected)

Determination (Unit)	LOI-371 %	LOI-425 %	LOI-650 %	LOI-1000 %
Assigned value	7.52	7.77	8.34	8.62
Number of Laboratories	5	8	11	17
Number of Analyses	30	44	54	80
Between Laboratory Standard Deviation	0.13	0.05	0.07	0.12
Within Laboratory Standard Deviation	0.06	0.03	0.04	0.10
Between Laboratory Relative Standard Deviation	1.68 %	0.69 %	0.86 %	1.35 %
Within Laboratory Relative Standard Deviation	0.79 %	0.35 %	0.45 %	1.19 %
Standard Uncertainty	0.051	0.018	0.021	0.019
Lower confidence limit	7.42	7.74	8.30	8.58
Upper confidence limit	7.62	7.81	8.38	8.65
Lower limit of tolerance	7.26	7.67	8.19	8.38
Upper limit of tolerance	7.77	7.88	8.48	8.85

Table 6 Fe-114 LOI Certified Values (non-hygroscopically corrected moisture)

Determination (Unit)	LOI-371 %	LOI-425 %	LOI-650 %	LOI-1000 %
Assigned value	7.76	8.06	8.58	8.85
Number of Laboratories	9	7	11	17
Number of Analyses	45	39	56	83
Between Laboratory Standard Deviation	0.15	0.10	0.08	0.13
Within Laboratory Standard Deviation	0.03	0.07	0.06	0.06
Between Laboratory Relative Standard Deviation	1.94 %	1.25 %	0.94 %	1.41 %
Within Laboratory Relative Standard Deviation	0.32 %	0.80 %	0.64 %	0.70 %
Standard Uncertainty	0.053	0.031	0.020	0.028
Lower confidence limit	7.66	8.00	8.54	8.80
Upper confidence limit	7.87	8.12	8.62	8.91
Lower limit of tolerance	7.46	7.86	8.42	8.60
Upper limit of tolerance	8.06	8.26	8.74	9.10

Origin of Material:

This material originates from the iron ore deposits of the Hamersley Group of the Hamersley Iron Province, central Pilbara Craton in Western Australia. The Hamersley Group rocks are early Proterozoic sediments, i.e. shale, siltstone, mudstone, chert, carbonate rocks and banded iron formations. Some minor volcanic units are also present in the Hamersley Group stratigraphy. There are two banded iron ore deposits within the Hamersley Group, the Brockman Iron Formation (approximately 2481±4 Ma) and the older Marra Mamba Iron Formation (approximately 2597±5 Ma). The Brockman Formation consists of cherts, mudstone, siltstone and banded iron formations; and is the main iron ore zone within the Hamersley Group. The iron formation consists of significant amounts of secondary enriched martite-goethite, which has an iron content of between 56 and 63 wt.%. The Marra Mamba Iron Formation is a surface enriched banded iron formation consisting of martite microplaty hematite ores, containing 60 to 68 wt.% iron.

1. Killick, M. F., H. M. Churchward, R. R. Anand (2003) Hamersley Iron Province, Western Australia. CRC LEME, Exploration and Mining.
2. Morris, R. C., M. Kneeshaw, (2011): Genesis modelling for the Hamersley BIF-hosted iron ores of Western Australia: A critical review. Australian Journal of earth sciences 58.

Minor Elements and Specific Gravity:

Informational values for minor elements are provided; results from Lazer Ablation ICP-MS and for Specific Gravity by Helium Pycnometer (results are for a single analysis at one laboratory):

SG (Helium Pycnometer) 4.32

Table 7 Fe-114 Minor Elements

As ppm	Ba ppm	Be ppm	Bi ppm	Cd ppm	Ce ppm	Co ppm	Cs ppm
18	54,5	0,8	0,16	-	27,1	7,5	0,11
Cu ppm	Dy ppm	Er ppm	Eu ppm	Ga ppm	Gd ppm	Ge ppm	Hf ppm
26	1,27	0,95	0,34	2,5	1,06	3,8	1,2
Ho ppm	In ppm	La ppm	Lu ppm	Mn ppm	Mo ppm	Nb ppm	Nd ppm
0,27	-	5,23	0,13	3090	1,2	1,81	5,53
Ni ppm	Pb ppm	Pr ppm	Rb ppm	Re ppm	Sb ppm	Sc ppm	Se ppm
22	65	1,61	0,4	-	1,2	6,5	-
Sm ppm	Sn ppm	Sr ppm	Ta ppm	Tb ppm	Te ppm	Th ppm	Tl ppm
1,36	1,2	5,2	0,63	0,17	0,4	1,93	-
Tm ppm	U ppm	V ppm	W ppm	Y ppm	Yb ppm	Zn ppm	Zr ppm
0,13	1,93	23,7	1,8	9	1,03	50	39,5

XRD Analysis:

A sample of Fe-114 was submitted for Quantitative X-Ray Diffraction Analysis after drying at 50°C.

XRDQUANT01 - Quantitative analysis, crystalline and amorphous content

Table 8 Fe-114 Mineral Abundance (XRD)

Mineral	% mass fraction
Hematite Fe_2O_3	23
Goethite $\text{FeO}(\text{OH})$	57
Magnetite Fe_3O_4	ND
Quartz SiO_2	2
Illite/Muscovite $\text{K}(\text{Al},\text{Mg},\text{Fe})_2(\text{Si},\text{Al})_4\text{O}_{10}(\text{OH})_2$	1
Amorphous content	17

Particle Size Distribution:

A sample of Fe-114 was submitted for particle size analysis by MICROTRAC S3500 Lazer Diffraction System.

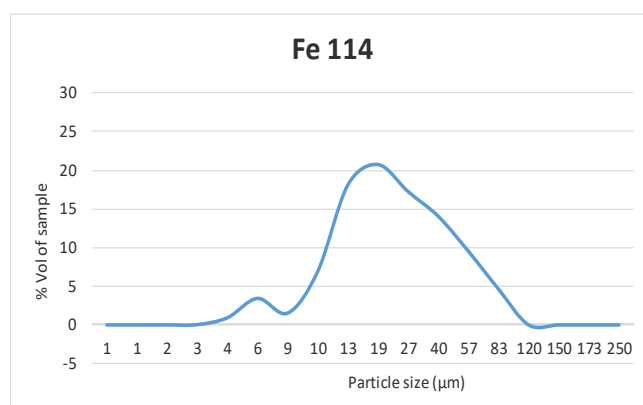


Figure 1 Fe-114 Particle Size Distribution

Spectral Analysis:

A sample of Fe-114 was submitted for analysis by TerraSpec 4 VNIR-SWIR.

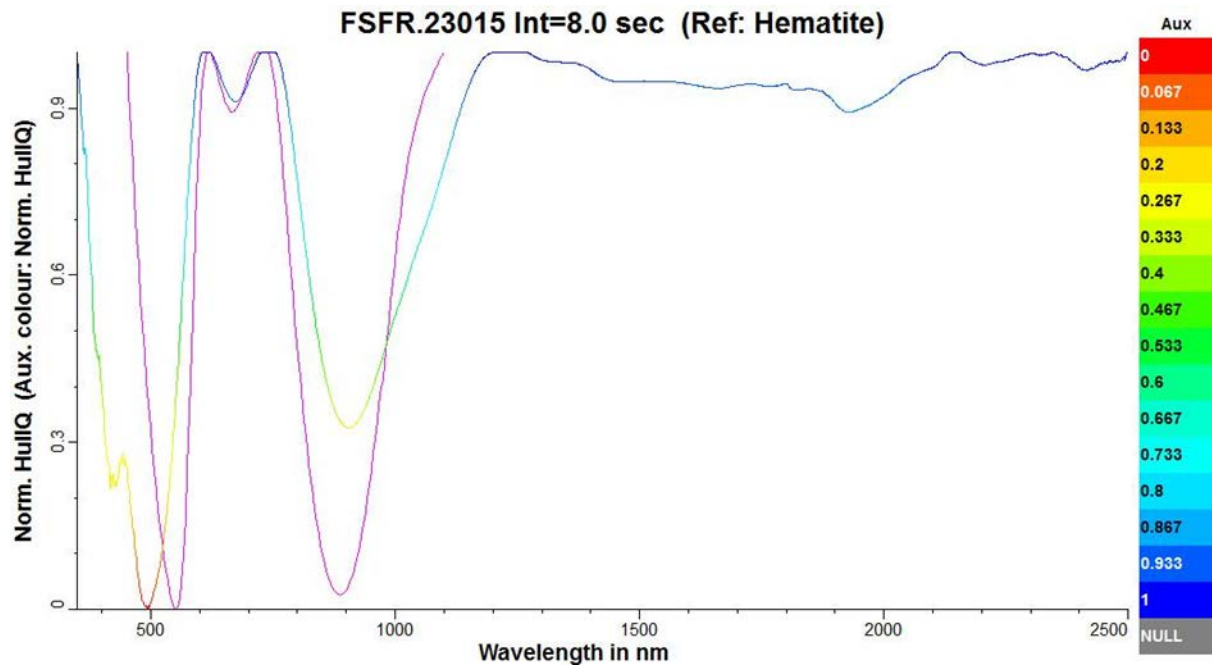


Figure 2 Fe-114 VNIR

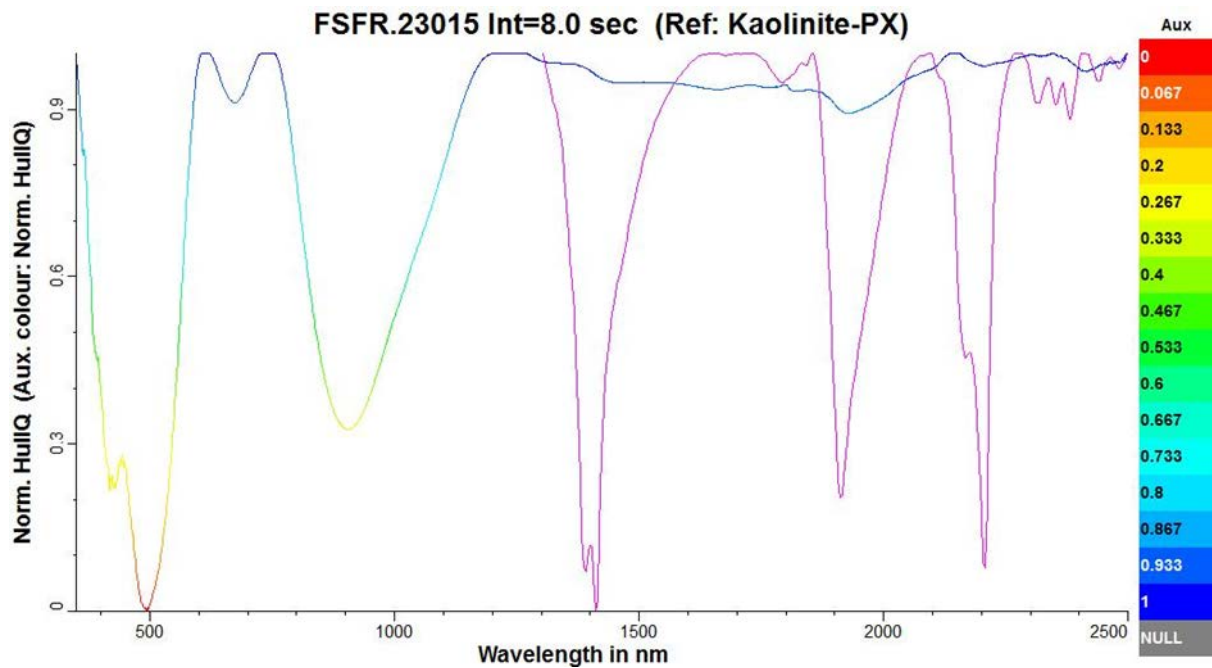


Figure 3 Fe-114 SWIR

FTIR Analysis:

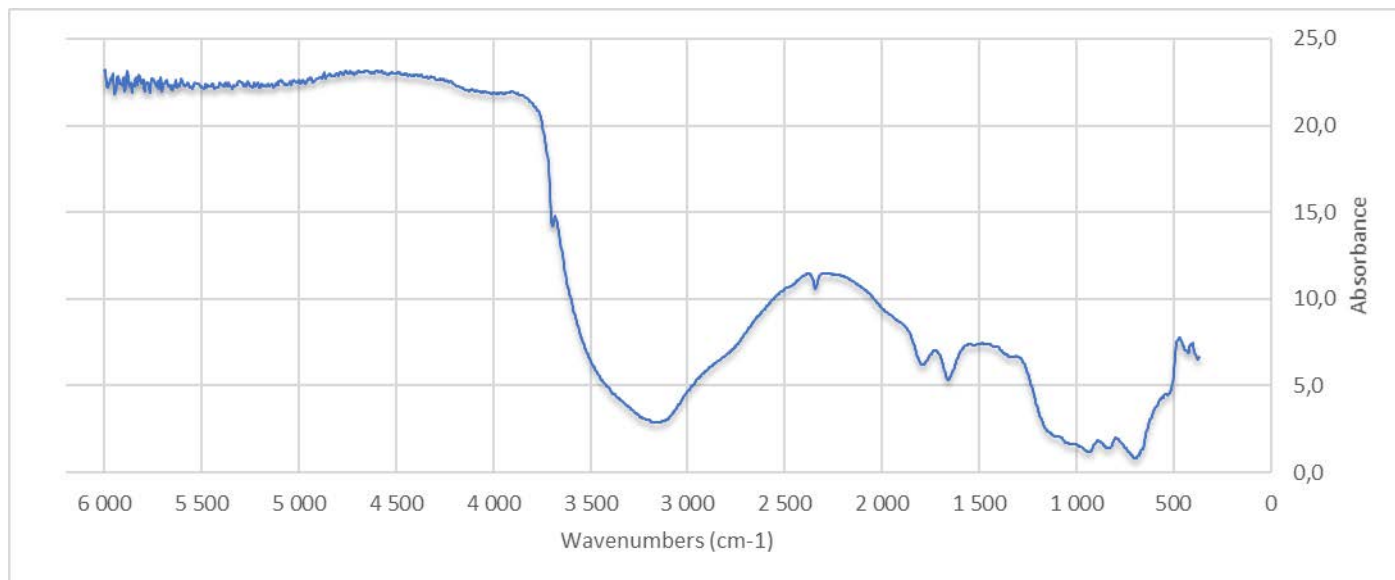


Figure 4 Fe-114 FTIR Spectra

Hygroscopic Moisture:

To illustrate the effect of hygroscopic moisture a sample of RRM Fe-114 underwent water vapour sorption testing.

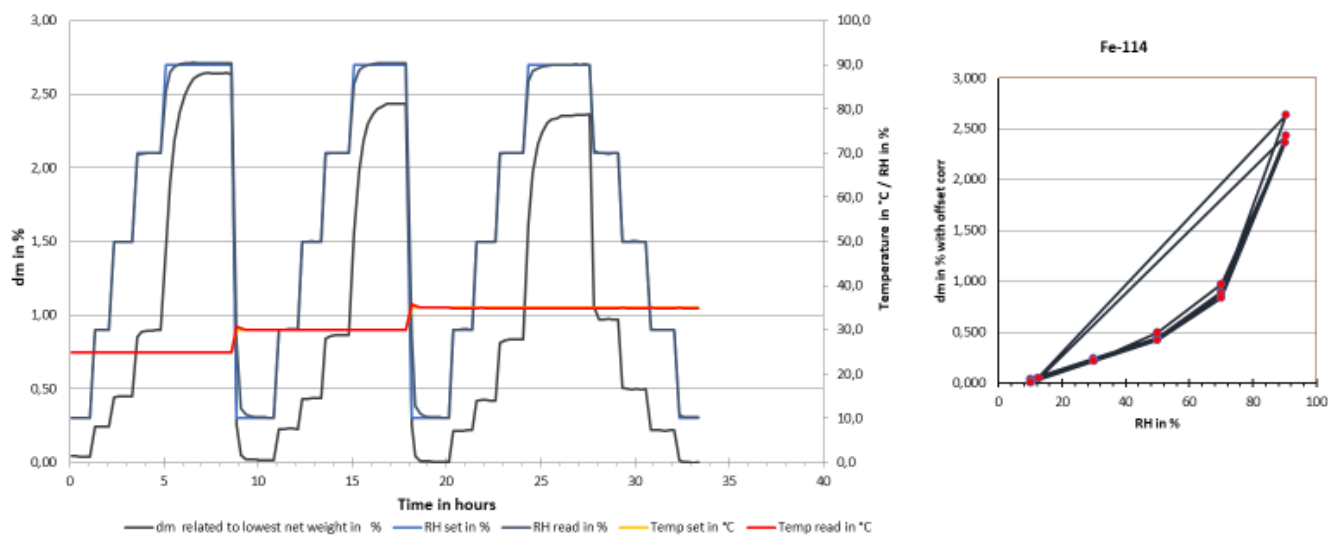


Figure 5 Fe-114 Sorption Testing

From the tests it is clear that hygroscopic moisture could have an impact on the quality of the chemical analysis results and appropriate storage and pre-treatment is highly recommended.

Participating Laboratories:

Table 9 Fe-114 Participating Laboratories

ALS Metallurgy, Australia
ALS Iron Ore, Australia
BHP Billiton Nelson Point Laboratory, Australia
BHP Billiton Newman Laboratory Line 1, Australia
BHP Billiton Newman Laboratory Line 2, Australia
BHP Billiton Newman Laboratory Line 3, Australia
BHP Billiton Newman Laboratory Line 4, Australia
Bureau Veritas, South Australia
Bureau Veritas, Western Australia
Citic Pacific Mining Laboratory, Australia
IMP Automation, Australia
Intertek Maddington, Australia
Intertek Robotic Laboratories Anderson Point, Australia
Intertek Robotic Laboratories Christmas Creek, Australia
Intertek Robotic Laboratories Cloudbreak, Australia
Intertek Robotic Laboratories FMG Port, Australia
Intertek Robotic Laboratories Mt Webber, Australia
Intertek Robotic Laboratories Roy Hill Port, Australia
Intertek Robotic Laboratories Solomon, Australia
Intertek, South Africa
Kumba Iron Ore Sishen Mine Laboratory, South Africa
MinAnalytical Laboratory Services, Australia
Rio Tinto Cape Lambert Port A Laboratory, Australia
Rio Tinto Cape Lambert Port B Laboratory, Australia
Rio Tinto Dampier Laboratory, Australia
Rio Tinto Marandoo Laboratory, Australia
Rio Tinto Paraburdoo Laboratory, Australia
Rio Tinto West Angelas Laboratory, Australia
Rio Tinto Yandicoogina Laboratory, Australia

Intended Use:

The Certified Reference Material RRM Fe-114 is intended to be used in analytical laboratories that analyse samples of similar grades (as presented in this certificate of analysis) and matrix.

- Monitoring of routine laboratory performance (both internal and external)
- Method development and method validation
- Instrument calibration

Instructions before Use:

Safety precautions (in-line with safe laboratory practices) for handling fine particulate matter are advised; such as the use of safety glasses, dust masks, gloves and laboratory coats.

Minimum Sample Size:

- Multi-element Iron Ore Suite – XRF fused disc analysis 0.7g
- LOI – 1g

Analysis Validity:

The property values for RRM Fe-114 Certified Reference Material remain valid provided that good laboratory practice is observed during handling and storage.

The material is in fine powder form and may be hygroscopic in nature.

Stability and storage:

This CRM is considered to have long-term stability under normal storage conditions. The CRM is available in different packaging and mass lots. When considering the lot size and said the ISO 16042, which recommends that the total contents of the container holding the CRM prior to its first use be divided in two portions that match the intended test portions.

Traceability:

The characterization of this material has been achieved by inter-laboratory study, each laboratory using an appropriate analytical method. These methods are calibrated against and analysed against high purity materials and appropriate Certified Reference Materials.

Additional Information:

Additional information on material RRM Fe-114 is available on request and includes the following:

- Data Tables
- Statistical Tables
- Particle Size Analysis Data
- Quantitative X-Ray Diffraction Analysis Report
- FTIR Spectral Data
- VNIR-SWIR Spectral Data
- Sorption Testing Data

Revisions:

This certificate is version 0.00 for the material RRM-Fe-114; and is intended to be a 'live document' intended to reflect progress in analytical chemistry. In that, any significant new data and information could be added at any time to ensure the currency and relevance of the certification. Any revisions to this Certificate of Analysis will be made available via the company website www.resourcereferencematerials.com

Legal Notice:

This material has been prepared and a comprehensive statistical evaluation conducted to assign the property values, with appropriate care and attention. The Purchaser, by receipt of this material and certificate, indemnifies and releases Resource Reference Materials (Pty) Ltd from and against any and all liability and costs arising from the use of this material and certificate and any actions taken thereupon.

Prepared and Supplied:

Certified Reference Material RRM Fe-114 has been prepared, certified and is supplied by;

Resource Reference Materials (Pty) Ltd
36 Michelson Road
Anderbolt, Boksburg
South Africa, 1459
E-Mail: info@resourcereferencematerials.com

Approved on behalf of Resource Reference Materials (Pty) Ltd:

Name
Position

C.S. van der Linde
Managing Director

29th May 2017



References:

ISO 17034:2016 General requirements for the competence of reference material producers

DRAFT ISO GUIDE 35:2016 Reference materials — Guidance for the characterization and the assessment of the homogeneity and stability of the material

ISO 5725-2:1994 Accuracy (trueness and precision) of measurement methods and results - Part 2: Basic method for the determination of repeatability and reproducibility of a standard measurement method

ISO 2596:2006 Iron ores — Determination of hygroscopic moisture in analytical samples — Gravimetric, Karl Fischer and mass-loss methods

ISO 11536:2015 Iron ores — Determination of loss on ignition — Gravimetric method

ISO 9516-1:2003 Iron ores — Determination of various elements by X-ray fluorescence spectrometry — Part 1: Comprehensive procedure

ISO 2597: 2006 part 3 Iron ores — Determination of hygroscopic moisture in analytical samples — Gravimetric, Karl Fischer and mass-loss methods

Copyright © Resource Reference Materials (Pty) Ltd